



- 1
- 2
- 3
- 4
- 5
- 6
- 7
- 8
- 9
- 10
- 11


FoodChain-Lab Installation

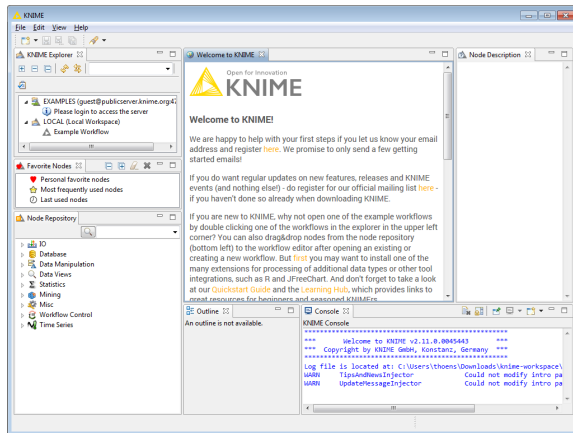


- The "KNIME Analytics Platform" must be installed before you can install FoodChain-Lab. Installation Instructions can be found at <https://tech.knime.org/installation-0>.
- When starting KNIME for the first time the interface looks like the image above.

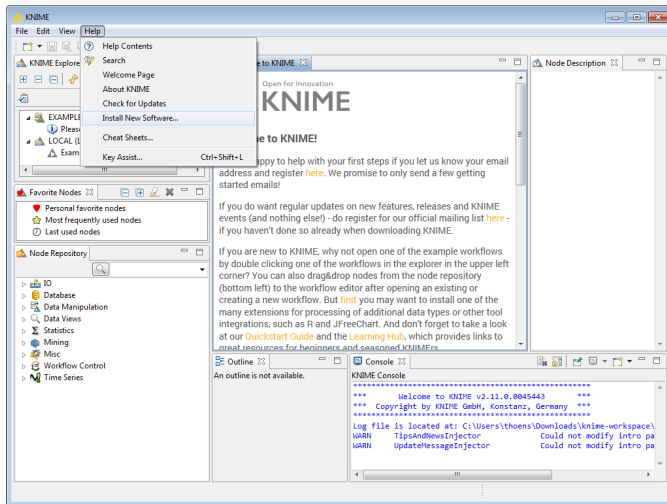


- Select **Help** > **Install New Software** in the menu bar.

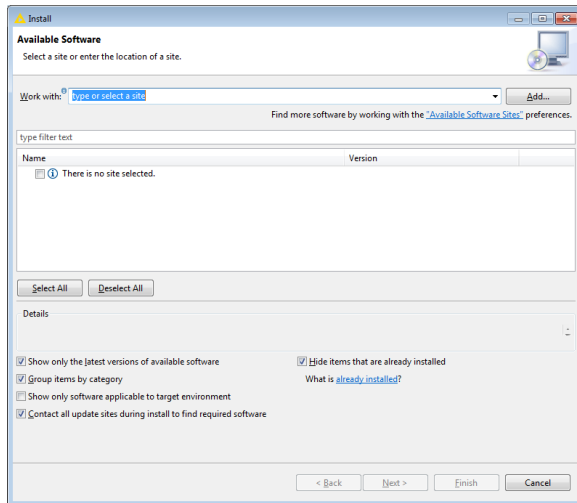


- Click **Add**, in the top-right corner.

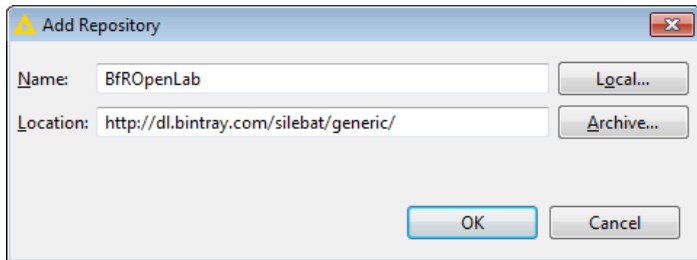


- In the Add Repository dialog that appears, enter "BfROpenLab" for the *Name* and the following URL for the *Location*:
`http://dl.bintray.com/silebat/generic/`
- Click **OK**.

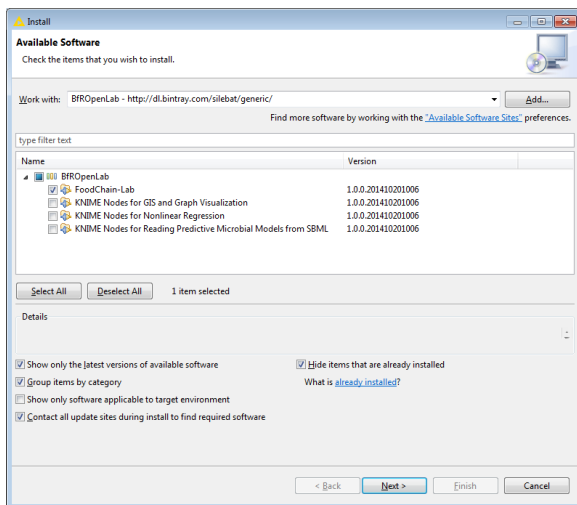


- In the Available Software dialog, expand the "BfROpenLab" entry and select the checkbox next to FoodChain-Lab.
- Click **Next**.

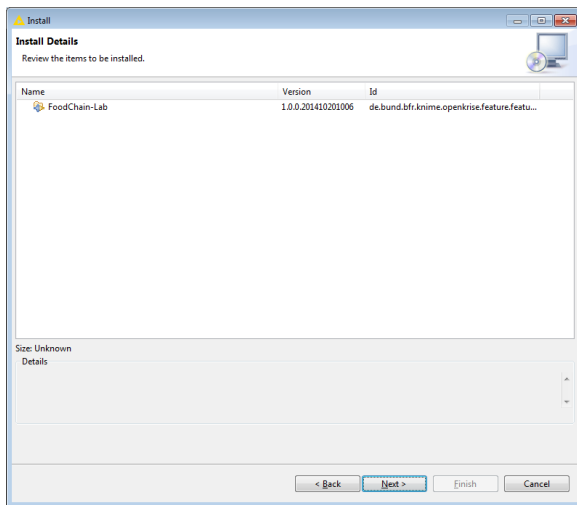


- In the next window, you'll see a list of the tools to be downloaded. Click **Next**.

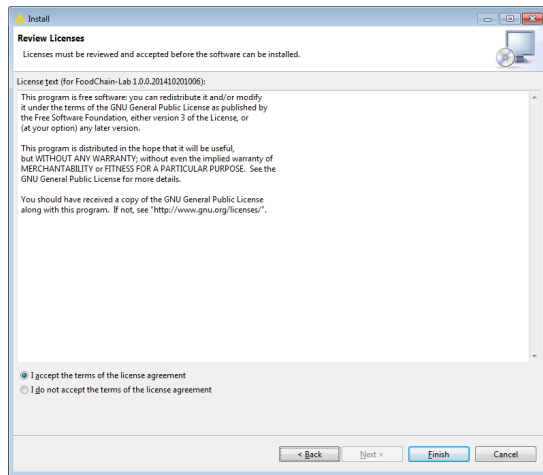


- Read and accept the license agreements, then click **Finish**.

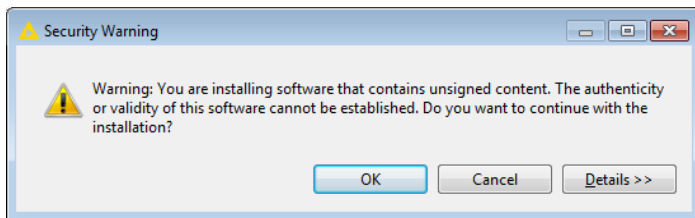


- If you get a security warning saying that the authenticity or validity of the software can't be established, click **OK**.

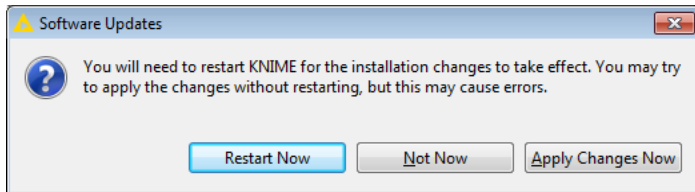

- When the installation completes, restart KNIME.

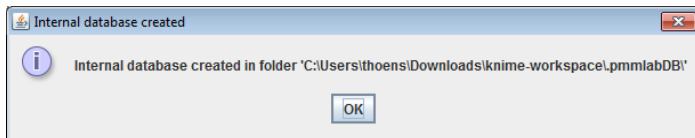

- When restarting you get a message saying that an internal database has been created. Click **OK**.

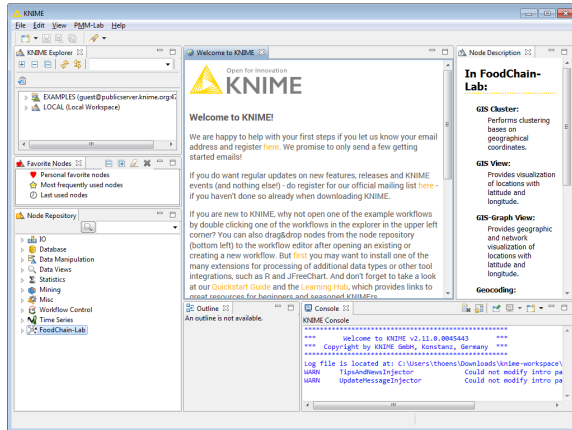

- When the KNIME interface has shown up, you should be able to see an item "FoodChain-Lab" in the **Node Repository** view in the bottom-left corner.